

## WEARABLE LIGHT THERAPY

### Patient Kit & Controller

CareWear delivers penetrating pulsed red and blue photonic energy, soothing warmth and ongoing relief from persistent pain and muscle tension. Your CareWear Patient Kit and Controller provides you with an excellent home-care maintenance option between in-clinic treatments and in conjunction with home exercise programs.

### Patient Kit Includes:

- User manual
- Quick start guide
- Single controller which can be used with all CareWear light patches
- Universal charger and USB 2.0 charger base
- Convenient zippered travel case with space for additional light patches



Practical,  
Easy-to-Use  
Patient Kit

## CAREWEAR COMPONENTS

### Reusable Light Patch



powerful LED light therapy in just 30 minutes.

Printed LEDs integrated into ultra-slim adhesive hydrogel patches, with up to 1500 LEDs per patch. Delivering

### Wireless Controller



is pre-programmed to deliver the correct dosage of light therapy to every CareWear light patch.

The controller acts as both power source and dosage control. Rechargeable and waterproof, the controller

### Replacement Light Patches



contains 10 patches which, with proper care, may be re-used up to 20 times each.

Replacement Light Patches are available in several body-contouring shapes and sizes. Each box

**CareWear**  
Wearable Therapeutics

© 2021 CareWear Corp. All rights reserved.  
1225 Financial Blvd. | Reno, NV 89509 | USA  
888-988-8797 | [www.carewear.net](http://www.carewear.net)

REORDER



Patents at: [www.carewear.net/patents](http://www.carewear.net/patents)

**CareWear**  
Wearable Therapeutics

## Wearable Light Therapy at Home

PATIENT GUIDE

Pain relief and recovery for:

HEAD/NECK

BACK

HEADACHE

TMJ

plus other areas as prescribed by your medical professional.

PhotoBioModulation and Soothing  
Warmth for Care Anywhere

# Pain Relief and Recovery Light Therapy for

HEAD/NECK

BACK

HEADACHE

TMJ

plus other areas as prescribed by your medical professional.

## Relieve Pain & Muscle Tension

The CareWear wearable light patch delivers penetrating pulsed red & blue photonic energy, soothing warmth & ongoing relief from persistent pain & muscle tension.

## Drug-Free Pain Relief

CareWear's Light Patches offer a drug-free alternative to Acetaminophen, NSAIDs and Opioids for controlling pain without the unwanted side effects.

### CareWear<sup>®</sup> Wearable Light vs. Pain Medication

	Decrease Pain	Decrease Inflammation	Safe for Long Term Use	No Side Effects
CareWear	✓	✓	✓	✓
NSAIDs	✓	✓		
Acetaminophen	✓			
Opioids	✓			



## On-the-Go Pain Relief

Accelerate recovery directly to the muscle or painful joint with CareWear's reusable adhesive light patches that provide ongoing relief from pain, while treating soft tissue injury.

## Improve Movement

- Reduce pain, joint stiffness and enhance circulation using CareWear Light Patches prior to exercising.
- Improve muscle endurance and joint stability.
- Get moving with drug-free Self Care and return to a healthier, more active lifestyle.

Use CareWear at home, between in-clinic treatments for accelerated recovery.

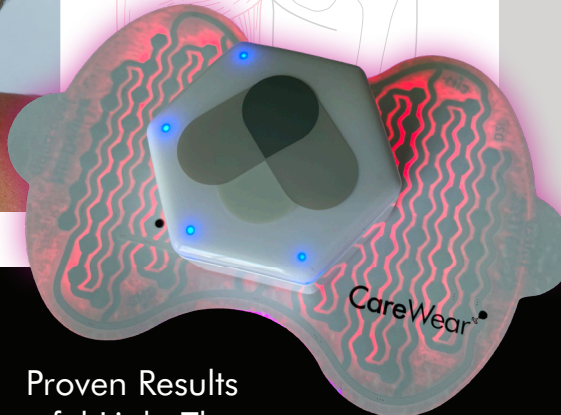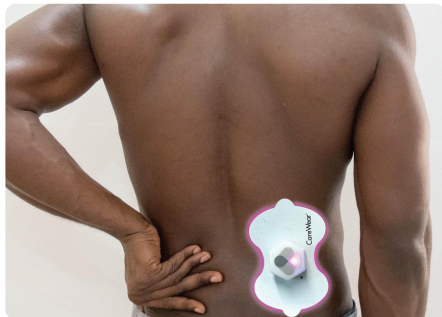

Proven Results  
Powerful Light Therapy

## FREQUENTLY ASKED QUESTIONS

### How many times can a light patch be worn?

The light patch use life will vary, depending on care and skin type, but should be replaced when it no longer adheres well to the skin. Correct skin preparation and proper patch storage can extend the use life of the light patch. The patches typically last for 20 uses.

### Can you make the light patch last longer?

Proper skin cleansing prior to patch placement can minimize build-up of sweat, oils, dead skin and hair on the patch's adhesive hydrogel. If the light patch has lost adhesion, rehydrate the hydrogel by moistening your finger with a few drops of water and gently rubbing the hydrogel to rehydrate. Proper patch handling to avoid air drying and UV light exposure can improve patch use life. After each treatment, be sure to place patch back on the clear liner and store in its resealable zip bag. Also avoid cutting, folding or pulling patch at non-tab points. Do not step on or provide excessive pressure to the patch socket.

### Can I use the CareWear light patch with other pain management therapies?

Yes. CareWear light patches are a non-drug therapy that can be used in combination with other therapies.

### Can I use the CareWear light patch with topical creams and gels?

Yes. However, topical creams/gels should NOT be applied to the treatment site immediately before applying the light patch as it may damage the hydrogel or produce excessive heat.

### Are CareWear light patches safe to use over tattooed skin?

The lower levels of fluence, irradiance and temperature change provided by the CareWear patches are insufficient to cause localized heating of red, blue or black tattoos.

### Can children use CareWear's light patches?

CareWear's light patches can be used on adults and children 13 years of age and older. Seek advice from a medical professional before applying to children under the age of 13.

CAREWEAR<sup>®</sup> TECHNOLOGY

Easy-to-Use, Unrestricted  
Delivery of Wearable Relief